

بوابة الخليج المبتكرة
Gulf Innovation Gate





Liner Bags

1.2.3 Liter



Canisters

1.2.3 Liter



Connection Pipe



T-Shaped Clamp



Gel Powder



Trolley





Wall Bracket




Regulator

Function Guide

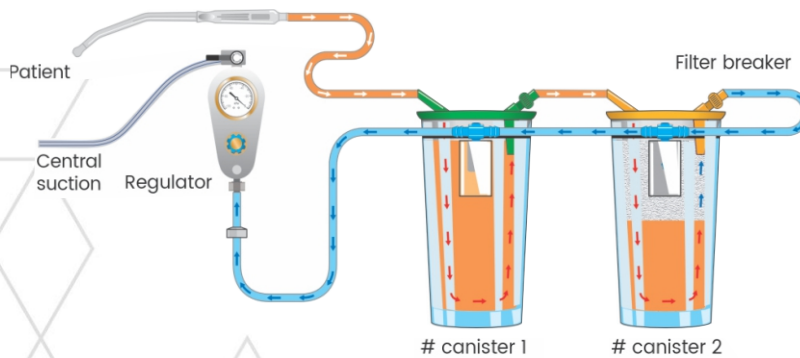

 The waste enters the primary bag and, as soon as the first bag is filled, passes through the connections to the next bag, which prevents the bag from overflowing because of a filter breaker.

 The suction system is connected to the canister or holder via the regulator, where the suction pressure is connected in a u-shape from the tank canister1 to the canister number 1. The suction pressure flows through the T-shaped connection between tanks 1 and 2, which is one of the features of Abadis Medical canisters which is a suitable suction pressure ratio. The filter breaker canister number 2 provides high reliability to prevent overflow.

Certified by:



بوابة الخليج المبتكرة
Gulf Innovation Gate



C.R No: 1437005

Contacts:

Address: Muscat, Sultanate of OMAN

Postal Code: 315 - P.O Box: 314

T e l : + 9 6 8 7 7 0 8 6 6 6 0

E-mail:Info@gulfinnovationgate.com

www.GULFINNOVATIONGATE.com